



V-Veil UP2™ Collector

Instructions for use

Intended purpose of V-Veil UP2™ Collector

V-Veil UP2™ Collector is a medical device used for the collection of secretions, various fluids, and substances from the body's natural cavities and from the skin.

Product overview

V-Veil UP2™ Collector is a disposable medical device, UV-C treated.

Materials

Polypropylene (PP)
Low-density polyethylene (LDPE)
100% Polyethylene/Polyester-bico

Patient information regarding self-sampling



Self-sampling is a method of specimen collection by the patient.

If you have questions about self-sampling, other specimen collection options, or for advice regarding preparation for collection, consult your healthcare provider, doctor, or pharmacist.

Sample collection

Important! Before collecting or self-collecting samples, follow the recommendations of the laboratory and/or healthcare professional regarding prior preparation, timing of collection, use during treatment or topical treatments, food intake, body cavity hygiene, smoking, local bleeding, etc.

Carefully read the instructions for use below. Use the device only once. The device must not be reused. For reliable results, follow the procedures carefully. If you have any questions about the procedure, contact your doctor or the laboratory.

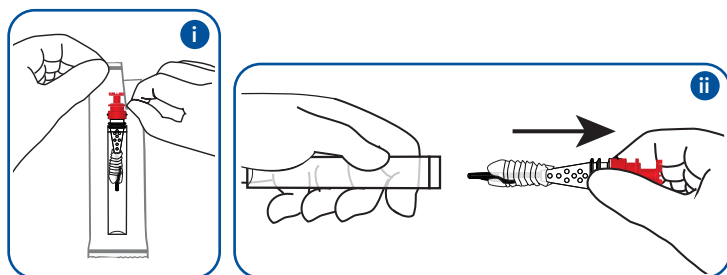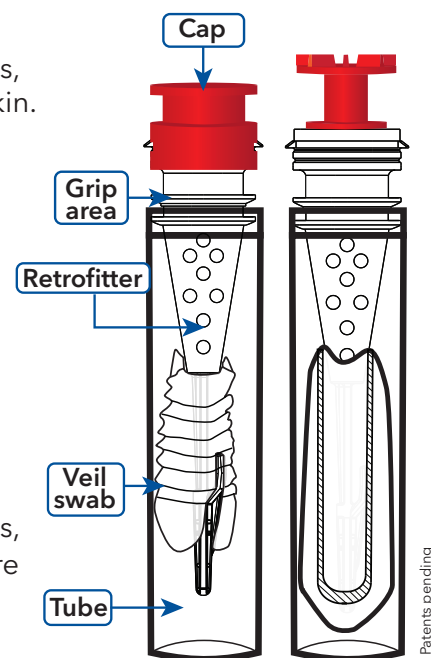


Open the packaging and remove the device (Diagram i). Do not touch the tip or veil swab - the soft white material! Remove the collection retrofitter from the tube, if provided with one (diagram ii), and place the tube on a flat, clean surface, following the instructions below to collect the specimen(s) requested.

Oral collection:


Follow the recommendations of the laboratory and/or health professional regarding the time of collection, food intake, mouth hygiene, smoking, etc.

Gently scrape the inside of the cheeks, under and on the tongue, rotating or flicking the swab to collect cells and/or saliva as shown in the adjacent pictures.

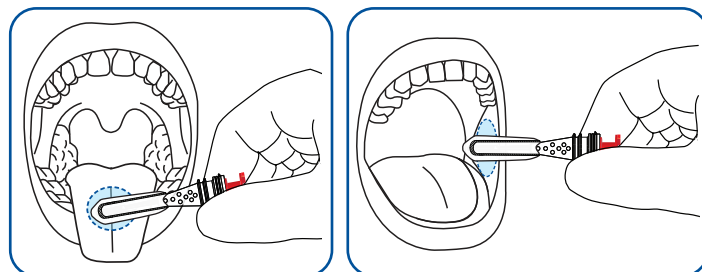

Anal collection:


The anal swab inspection method is used to check for the presence of parasites or other pathogens. Follow the recommendations of the laboratory and/or healthcare specialist. Collect the specimen by simply swabbing or rubbing the outer part of the anus with the soft white part of the device, or gently insert the veil swab into the anus for a distance of 3-5 cm, then rotate gently to remove it.

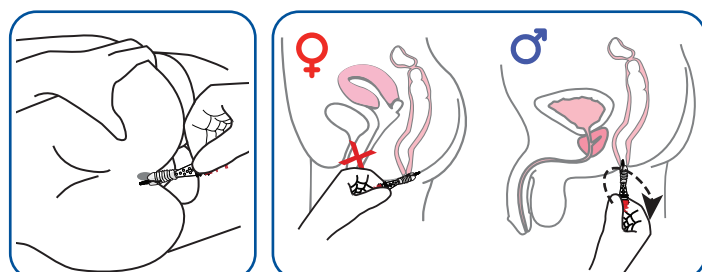

For **collection from the skin** rub the veil swab over the targeted area (groin, perianal, axial, etc.), applying moderate pressure, to collect cells and secretions.

Penile meatus specimen collection:

Hold the penis with your free hand (the hand without the device). Using the other hand (with the retrofitter), gently insert the first 5 millimeters of the thin end into the external urethral opening - the opening of the penis through which you urinate - and rotate gently clockwise once, only if recommended. Ensure you rotate completely to obtain the best specimen.

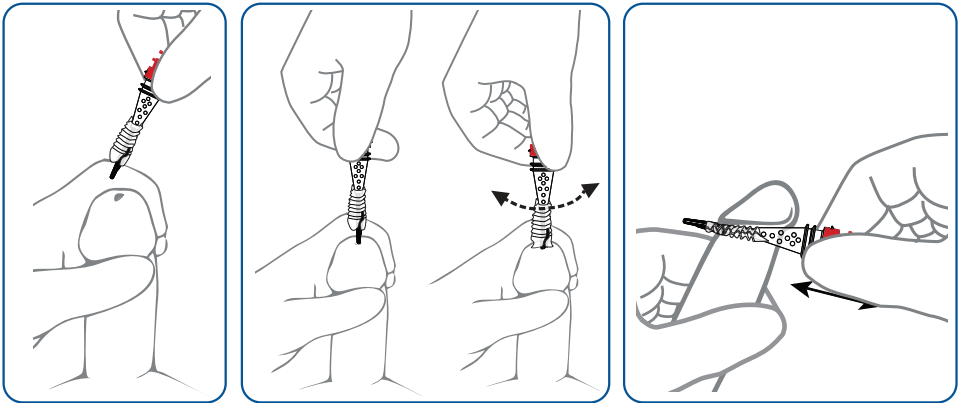


It is not necessary to insert deeply into the opening of the penis, nor to rotate with too much force so as to cause injury. After collection from the urethra (if recommended), wipe the surface under the glans, the balanopreputial sulcus area, using the part provided with the veil - the soft white material.

Insert the collected sample into the tube (if provided) and press to close the cap (Diagram iii). Return the medical device with the collected sample as instructed by your doctor or laboratory.

Storage

Store in a clean, dry place away from heat sources, direct sunlight and water. After collection, it is recommended to close the cap and store the specimen as advised by your doctor or laboratory, in a clean, dry place, away from heat sources and direct sunlight.

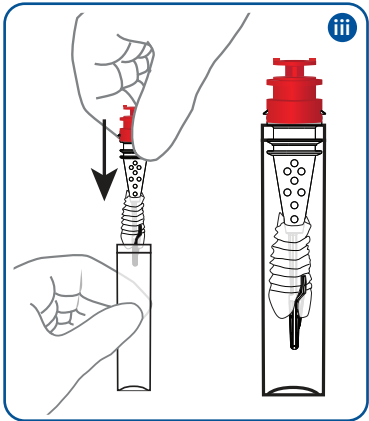

Incident reporting

IMPORTANT! Only the instructions for use may be used as reference for legal information, recommendations and restrictions. Do not use the product if the packaging is damaged and if the expiration date has passed. If, despite the measures we take during the manufacturing process, you believe you have received a defective product, please do not use it and contact the manufacturer or distributor, clearly stating the problem encountered and the batch number specified on the packaging. Device-related incidents can be reported to the manufacturer at client@vvup.org and to the competent authority of the Member State where you are located. The V-Veil UP2™ Collector and requirements of the European Union.

CE	European Conformity
MD	Medical Device
	UV-C treated
	Do not reuse
	Attention
	Read the notice
	Do not use if package is damaged
	Keep dry
	Keep away from direct sunlight

Authorized importer (if applicable): see information on the packaging
 Authorized distributor (if applicable): see information on the packaging

V-VEIL-UP PRODUCTION SRL
Ciocanesti 280B Calinesti
117197 Arges ROMANIA
client@vvup.org
TEL : 0040 348 420 620