



V-Veil UP2™ Anal

Instructions for use

Intended purpose of V-Veil UP2™ Anal

V-Veil UP2™ Anal is a medical device used for the collection of secretions, various fluids, and substances from the anal cavity.

Product overview

V-Veil UP2™ Anal is a disposable medical device, UV-C treated.

Materials

Polypropylene (PP)

Low-density polyethylene (LDPE)

100% Polyethylene/Polyester-bico

Recommendations

To avoid sample contamination, it is recommended:

- Avoid receptive sexual contact, at least 24-48 hours beforehand.
- Do not use topical products: creams, ointments, gels, lubricants or suppositories in the anal area.
- Do not perform douches (enemas): avoid rectal showers or enemas in the last 24 hours.
- For women during their menstrual period: menstruation does not directly affect the anal test, but it is mandatory to inform the doctor if you use an internal tampon during the procedure.

Patient information regarding self-sampling

Self-sampling is a method of specimen collection by the patient.

If you have questions about self-sampling, other specimen collection options, or for advice regarding preparation for collection, consult your healthcare provider, doctor, or pharmacist.



For the self-sampling of anal specimens in minor patients, it is recommended to consult local, national, and federal regulations.


Anal specimen collection

The anal cavity specimen collection device can be used by both women and men. Wash your hands before starting.

Open the device packaging (diagram 1). Remove the retrofitter from the tube (if provided with one) as shown in diagram 2 and place the tube on a clean, flat surface before performing the collection. Do not touch the tip or the veil swab - the soft white material - and do not place the retrofitter on any surface! If the tip or the veil swab - the soft white material - is touched, placed down, or dropped, use a new collection device. Holding the grip area/end of the retrofitter, carefully insert the veil swab part - the soft white material (the veil) - into the rectum approximately 3-5 cm past the anal verge (the outer part of the anus), according to diagram 3 for women and diagram 4 for men.

After insertion, gently rotate clockwise 2 times (diagram 5) and counterclockwise another 2 times (diagram 6). The entire collection procedure will take from 5 to 10 seconds.

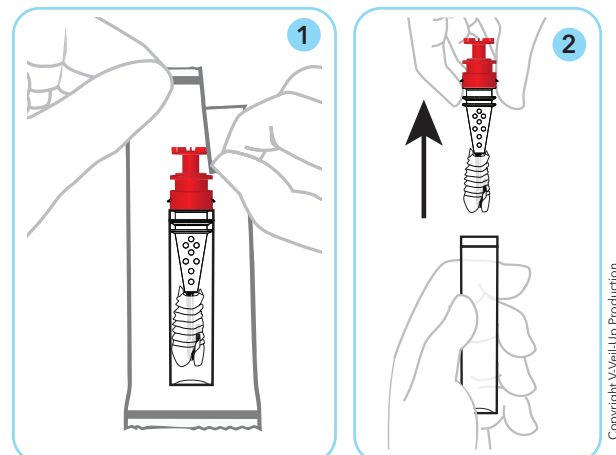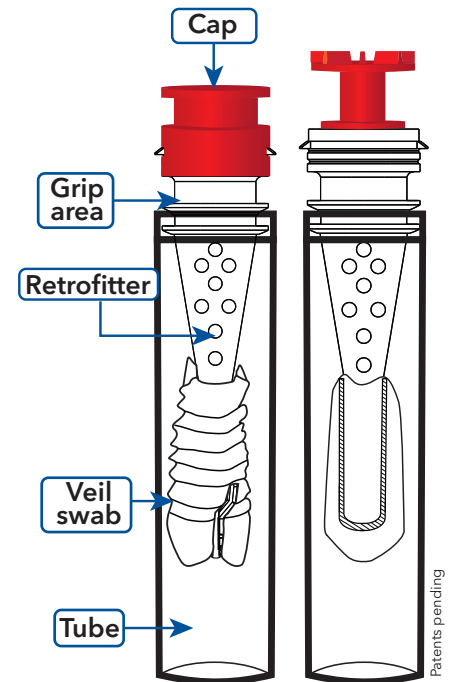

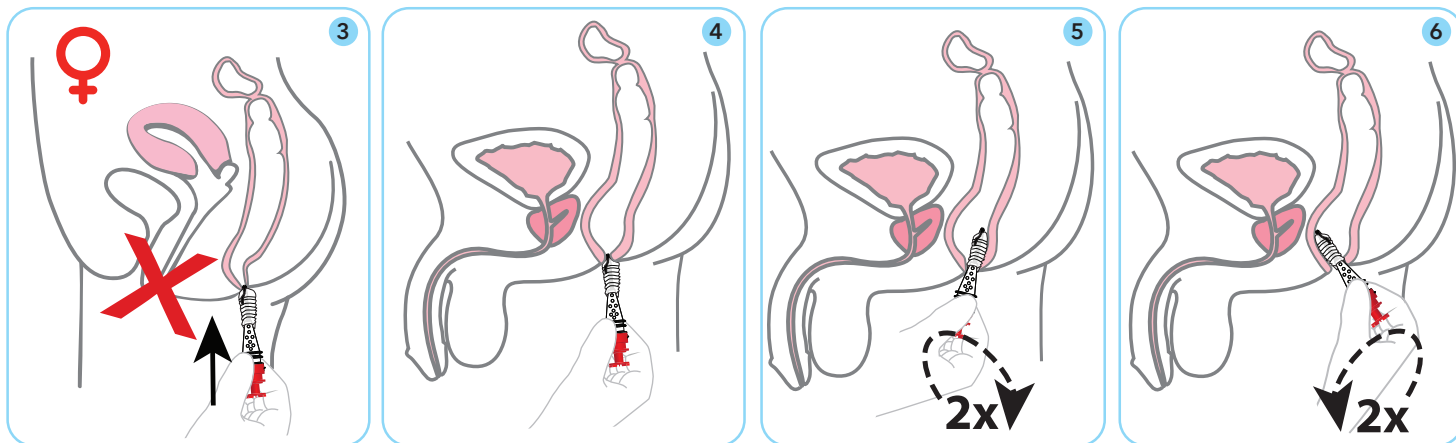


Withdraw the retrofitter gently, without touching the skin (diagram 7). After collection, insert the specimen into the tube (if provided) and press to close the cap, according to diagram 8. Return the medical device with the collected specimen (diagram 9) according to your doctor's or laboratory's instructions.

Storage

Store in a clean, dry place away from heat sources, direct sunlight and water.

After collection, it is recommended to close the cap and store the specimen as advised by your doctor or laboratory, in a clean, dry place, away from heat sources and direct sunlight.

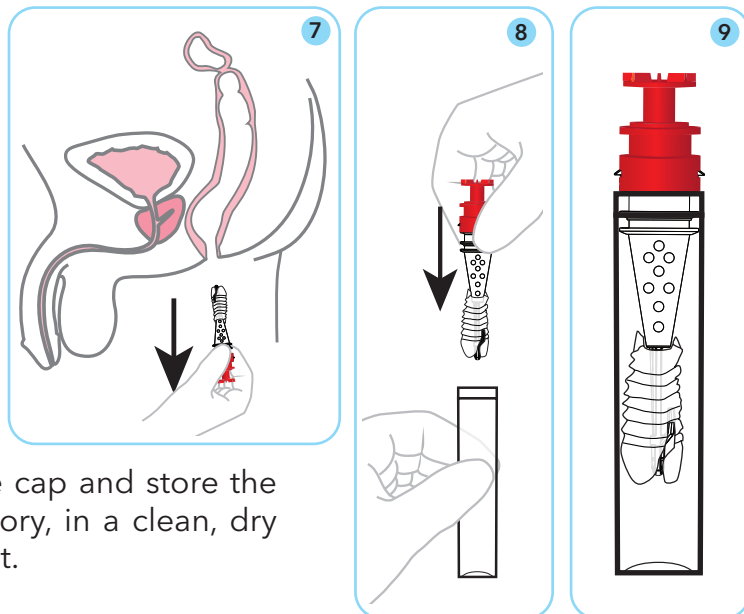

Incident reporting

IMPORTANT! Only the instructions for use may be used as reference for legal information, recommendations and restrictions. Do not use the product if the packaging is damaged and if the expiration date has passed. If, despite the measures we take during the manufacturing process, you believe you have received a defective product, please do not use it and contact the manufacturer or distributor, clearly stating the problem encountered and the batch number specified on the packaging. Device-related incidents can be reported to the manufacturer at client@vvup.org and to the competent authority of the Member State where you are located. The V-Veil UP2™ Anal and requirements of the European Union.

	European Conformity
	Medical Device
	UV-C treated
	Do not reuse
	Attention
	Read the notice
	Do not use if package is damaged
	Keep dry
	Keep away from direct sunlight

The instructions for use are available for the following codes:

ANR60V14, ANR93V24, ANR60V13F, ANR93V23F, ANR60V134, ANR60V, ANR93V.



Authorized importer (if applicable): see information on the packaging



Authorized distributor (if applicable): see information on the packaging