


## Simultaneous sampling of cells from the endocervix and ectocervix

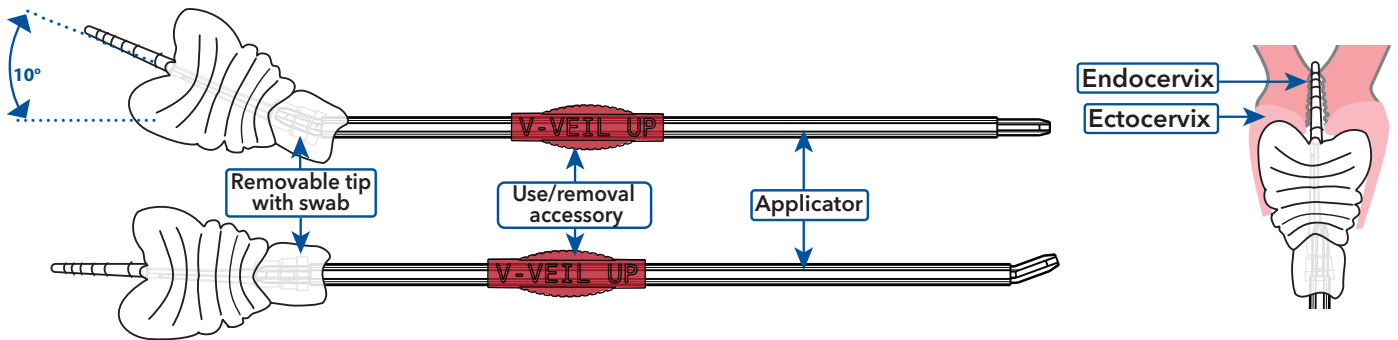

### Intended purpose of V-Veil UP2™ Cervical Brush

V-Veil UP2™ Cervical Brush is a medical device used for the collection of cervical cells, secretions, various fluids and substances from the cervix and vaginal cavity.

### Product overview

Disposable medical device, UV-C treated, with removable tip. For use by qualified medical personnel only. The medical device V-Veil UP2™ Cervical Brush is made of Polypropylene (PP and LDPE) and 100% Polyethylene/Polyester-bico. The device can be provided with a use/removal accessory.

### Contraindications

The device should not be used during pregnancy.

Do not use for people allergic to plastics.


Do not use lubricants or other products containing substances such as carbomers.

### Sample collection



**Important!** Before collecting or self-collecting samples, follow the recommendations of the laboratory and/or healthcare professional regarding prior preparation, timing of collection, use during treatment or topical treatments, food intake, body cavity hygiene, smoking, local bleeding, etc.

Carefully read the instructions for use below. Use the device only once. The device must not be reused. For reliable results, follow the procedures carefully. If you have any questions about the procedure, contact your doctor or the laboratory.

Open the packaging and remove the device (diagram 1). Depending on the position of the cervix, you may choose to use the 10° angled end device shown in Diagram 2 or the non-angled end device shown in Diagram 3. Insert the tip of the brush into the endocervical canal deep enough to allow the veil to fully contact the ectocervix. Gently press and rotate the brush at least 2 times clockwise and 2 times counterclockwise forming at least 4 complete 360° rotations, then retract it.

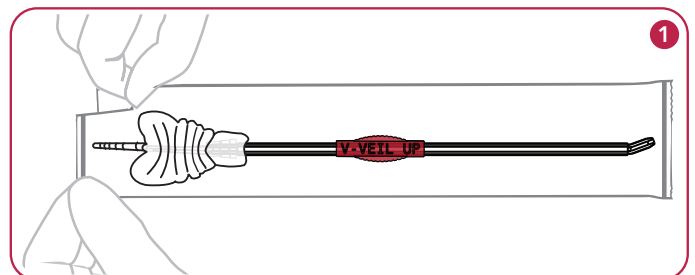

**WARNING:** If the device is fitted with the use/removal accessory, be careful not to insert the accessory inside the vaginal cavity!

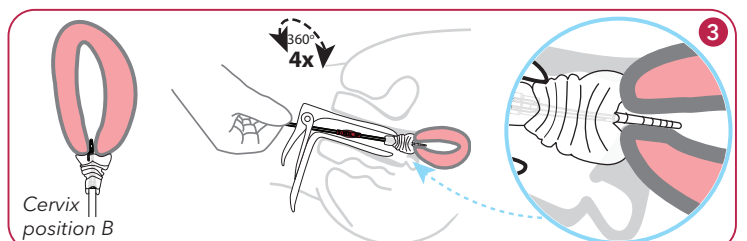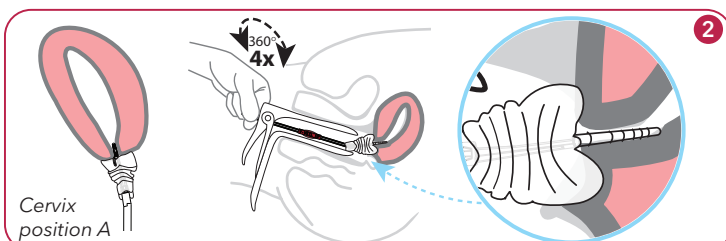


Preparation

Option A.


Place the device with the sample collected as quickly as possible in the container required for sending to the laboratory. With the device completely immersed in the preservation solution or biological sample fixative, rotate it vigorously at least 2 to 6 times in each direction to release the material, as shown in Diagram 4. Visually inspect the V-Veil UP2™ Cervical Brush to make sure there is no material left attached, and if there is, repeat the procedure in Diagram 4.

Depending on the manufacturer’s recommendations of the vial and container used, you can discard the collection device after insertion (diagram 5a) or leave it inside (diagram 5b). The device does not float (diagram 6). Store and return the sample according to your laboratory’s instructions.

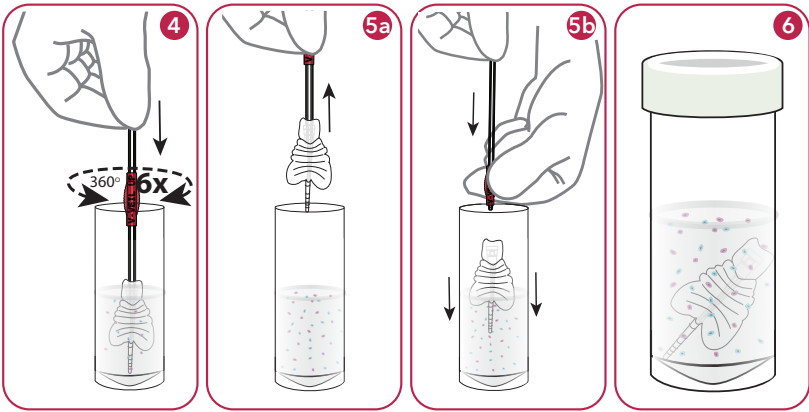


**Option B.** Transfer the cellular material onto the slide using painting-like motion, first on one side and then on the other side of the device (diagram 7). Immediately apply fixative to the slide (diagram 8).

ATTENTION! Collection and transportation container not included.

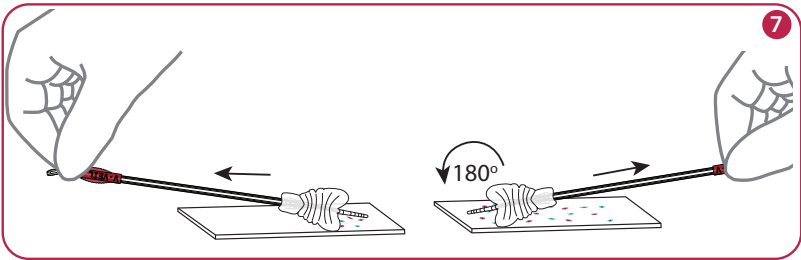

Storage

Store in a clean, dry place away from heat sources, direct sunlight and water. After collection, it is recommended to store the sample according to the instructions of your laboratory, preservation solution or biological sample fixative. If the collected sample is stored incorrectly, its effectiveness may be compromised.

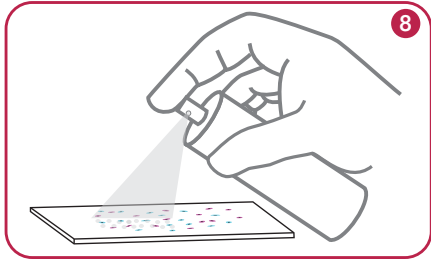

Incident reporting

IMPORTANT! Only the instructions for use may be used as reference for legal information, recommendations and restrictions. Do not use the product if the packaging is damaged and if the expiration date has passed. If, despite the measures we take during the manufacturing process, you believe you have received a defective product, please do not use it and contact the manufacturer or distributor, clearly stating the problem encountered and the batch number specified on the packaging. Device-related incidents can be reported to the manufacturer at [client@vvup.org](mailto:client@vvup.org) and to the competent authority of the Member State where you are located. The V-Veil UP2™ Cervical Brush is manufactured in compliance with the applicable provisions and requirements of the European Union.

BRS198V7635C

BRS202V7635

|  |                                  |
|--|----------------------------------|
|  | European Conformity              |
|  | Medical Device                   |
|  | UV-C treated                     |
|  | Do not reuse                     |
|  | Attention                        |
|  | Read the notice                  |
|  | Do not use if package is damaged |
|  | Keep dry                         |
|  | Keep away from direct sunlight   |

Patents pending

Authorized importer (if applicable): see information on the packaging  
 Authorized distributor (if applicable): see information on the packaging

REF IFU: BRS-EN V02 2026.10

V-VEIL-UP PRODUCTION SRL

Ciocanesti 280B Calinesti  
117197 Arges ROMANIA  
[client@vvup.org](mailto:client@vvup.org)  
TEL : 0040 348 420 620