

Intended purpose of V-Veil UP2™ Oral

V-Veil UP2™ Oral is a medical device used for the collection of oropharyngeal and oral secretions by deep expectoration and/or spitting.


Product overview

Disposable medical device, UV-C treated. The medical device V-Veil UP2™ Oral is made of Polypropylene (PP) and 100% Polyethylene/Polyester-bico. The device can be fitted with fasteners so that it can be used with various tubes.

Sample collection

Follow the recommendations of the laboratory and/or health professional in terms of sample preparation, timing, food ingestion, mouth hygiene, smoking, etc.

Open the packaging and remove the device (diagram 1). Grasp the veil holder (the oval piece of hard material) by the two protruding parts on the outside, holding it by the grip area (Diagram 2). Cough vigorously from the inside of the chest for 3-6 minutes on the veil - the soft white material as shown in Diagram 3. If your doctor or specialist has recommended that you spit out your sputum to identify a possible bacterial, viral, or fungal infection, you can do so on a veil. If your doctor or specialist has recommended a successive use, strictly follow the instructions of qualified medical personnel.

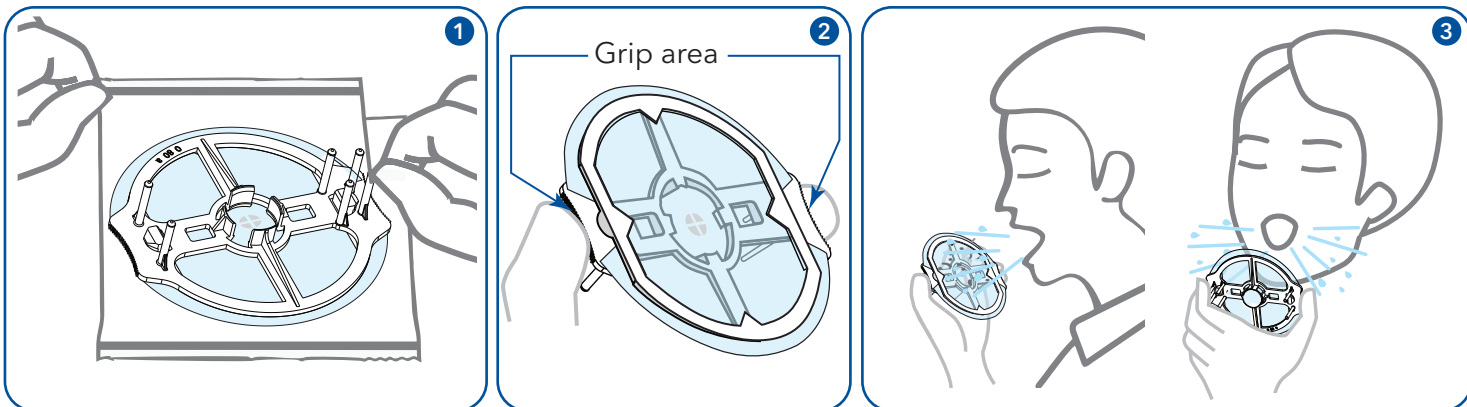

ATTENTION! Collection and transportation container not included. Use the container recommended or provided by your healthcare professional.


Storage

Store in a clean, dry place away from heat sources, direct sunlight and water. After collection, it is advisable to store the sample according to the instructions of qualified medical personnel and according to the recommendations of the preservation solution or biological sample fixative, if used. If the collected sample is stored incorrectly, its effectiveness may be compromised.

Incident reporting

IMPORTANT! Only the instructions for use may be used as reference for legal information, recommendations and restrictions. Do not use the product if the packaging is damaged and if the expiration date has passed. If, despite the measures we take during the manufacturing process, you believe you have received a defective product, please do not use it and contact the manufacturer or distributor, clearly stating the problem encountered and the batch number specified on the packaging. Device-related incidents can be reported to the manufacturer at client@vvup.org and to the competent authority of the Member State where you are located. The V-Veil UP2™ Oral is manufactured in compliance with the applicable provisions and requirements of the European Union.

 Authorized importer (if applicable): see information on the packaging
 Authorized distributor (if applicable): see information on the packaging

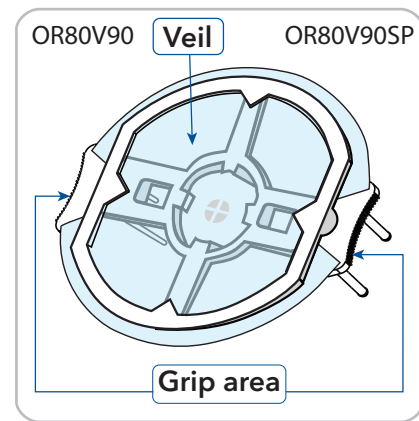

Patents pending

	European Conformity
	Medical Device
	UV-C treated
	Do not reuse
	Attention
	Read the notice
	Do not use if package is damaged
	Keep dry
	Keep away from direct sunlight