


V-Veil UP2™ Tiny Swab

Instructions for use

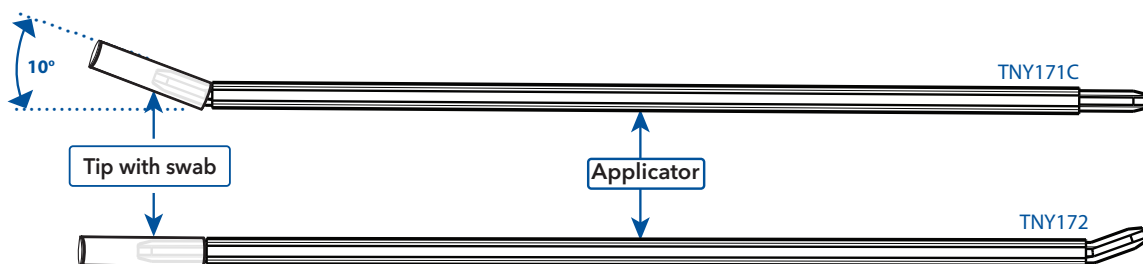

Intended purpose of V-Veil UP2™ Tiny Swab

V-Veil UP2™ Tiny Swab is a medical device used for the collection of secretions, various fluids and substances from natural body cavities and skin.

Product overview

Disposable medical device, UV-C treated. The medical device V-Veil UP2™ Tiny Swab is made of Polypropylene and 100% Polyethylene/Polyester-bico.

Contraindications

Do not use for people allergic to plastics.

Patient information on self-collection

Self-collection is a way of collecting samples by the patient.

If you have questions about self-collection, other sample collection options, and need advice on how to prepare for collection, consult your healthcare provider, doctor or pharmacist.

Sample collection


Important! Before collecting or self-collecting samples, follow the recommendations of the laboratory and/or healthcare professional regarding prior preparation, timing of collection, use during treatment or topical treatments, food intake, body cavity hygiene, smoking, local bleeding, etc.

Carefully read the instructions for use below. Use the device only once. The device must not be reused. For reliable results, follow the procedures carefully. If you have any questions about the procedure, contact your doctor or the laboratory.

Open the packaging and remove the device (Diagram i). Do not touch the tip or veil swab - the soft white material and do not place the device on any surface!

Nasal collection


Carefully insert the swab into one nostril, advancing it along the lower nasal passage until you feel a slight resistance. Roll the swab over the mucosa for a few seconds to collect a sufficient amount of secretion. Use the same swab to repeat the procedure in the other nostril.

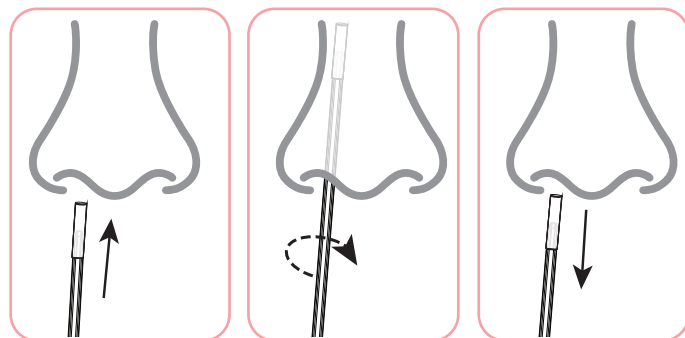

Oral collection:

Follow the recommendations of the laboratory and/or health professional regarding the time of collection, food intake, mouth hygiene, smoking, etc.

Gently scrape the inside of the cheeks, under and on the tongue, rotating or flicking the swab to collect cells and/or saliva as shown in the adjacent pictures.

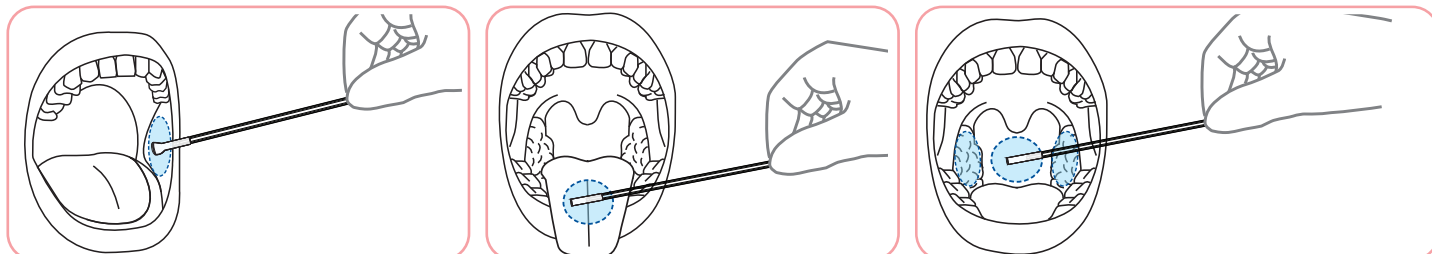


Vaginal collection

The device should not be used during pregnancy.

Do not use lubricants or other products containing substances such as carbomers.

The sample will be collected before the pelvic examination, using a speculum or manually.

Gently insert the veil swab into the vagina as shown in Diagram a. If the veil does not slip, turn it gently as you push. If it is still difficult, stop.

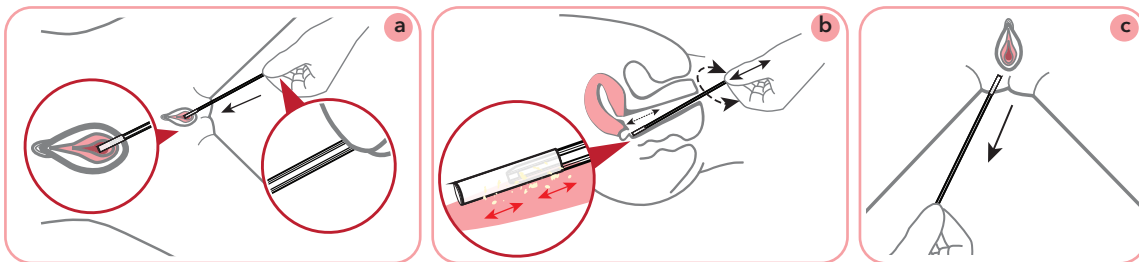


Rotate or move the veil swab back and forth for 10 to 15 seconds, as shown in Diagram b (using the use/removal accessory, if provided). Remove the swab without touching the skin outside the vagina, as shown in Diagram c.

Anal collection

The anal swab inspection method is used to check for the presence of parasites or other pathogens.

The patient is placed in a knee-chest position or lateral decubitus position. Gently insert a veil swab 3-5 cm into the anus and rotates 4-6 times before removing it.

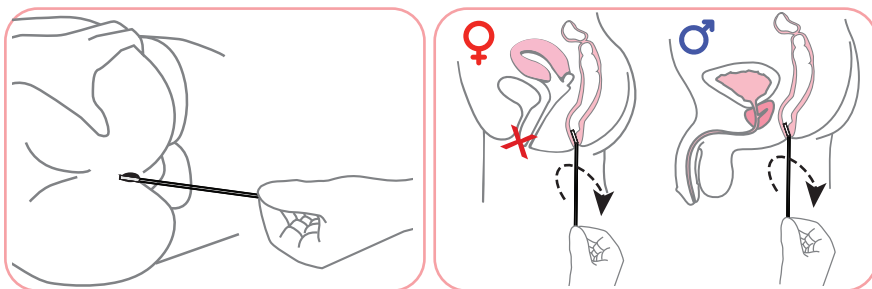


For **collection from the skin** rub the veil swab over the targeted area (groin, perianal, axial, etc.), applying moderate pressure, to collect cells and secretions.

Option A After collection, insert the device with the collected sample into the container required for sending to the laboratory. Using scissors, cut the applicator, approximately 1 cm from the end with the veil, leaving the part with the impregnated veil to fall freely into the collection container.

WARNING! Collection and transportation container not included.

Option B Transfer the cellular material onto the slide using painting - like motion, first on one side and then on the other side of the device. Immediately apply fixative to the slide.

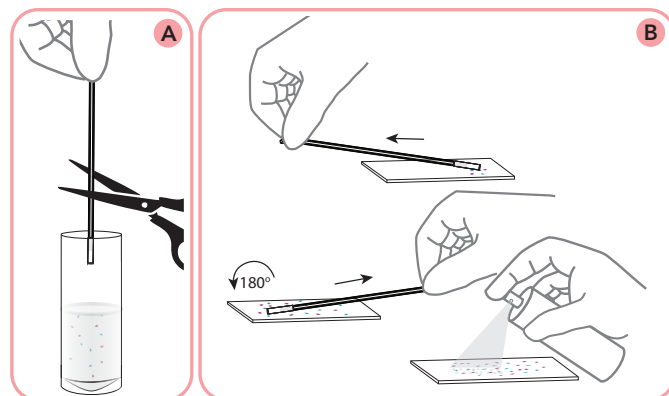

Storage

Store in a clean, dry place away from heat sources, direct sunlight and water. After collection, it is recommended to close the lid of the collection container and store the sample as recommended by the laboratory or health professional.

Incident reporting

IMPORTANT! Only the instructions for use may be used as reference for legal information, recommendations and restrictions. Do not use the product if the packaging is damaged and if the expiration date has passed. If, despite the measures we take during the manufacturing process, you believe you have received a defective product, please do not use it and contact the manufacturer or distributor, clearly stating the problem encountered and the batch number specified on the packaging. Device-related incidents can be reported to the manufacturer at client@vvup.org and to the competent authority of the Member State where you are located. The V-Veil UP2™ Tiny Swab is manufactured in compliance with the applicable provisions and requirements of the European Union.

CE	European Conformity
MD	Medical Device
UV-C TREATED	UV-C treated
Do not reuse	Do not reuse
Attention	Attention
Read the notice	Read the notice
Do not use if package is damaged	Do not use if package is damaged
Keep dry	Keep dry
Keep away from direct sunlight	Keep away from direct sunlight

Patents pending